



CERTIFICATION by CUSTODIAN of RECORDS
Tucson Police Department Crime Laboratory

PK
7.10.2007

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002265 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 7-5-07 Cal 1 0.101 AC Cal 2 0.099 AC

Between the dates of 6-7-07 and 7-5-07 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.
The instrument was returned to service on N/A.


Clare Conley

Clare L. Conley

7-9-07

ADPS QAS Permit # 33314

SUBSCRIBED AND SWORN to before me this 9, day of July, 2007.

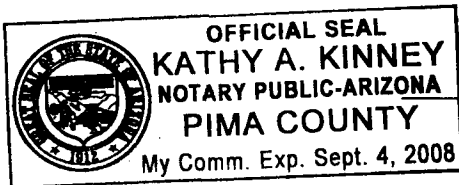

Kathy A Kinney
Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

TS4

QA SPECIALIST Clare Conley AGENCY Tucson Police Department
DATE 7-5-07 TIME 1250
INTOXILYZER SERIAL # 80-002265 LOCATION ODD

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
(✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
(✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
(✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
(✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 099 AC

Instrument is operating properly and accurately. YES ✓ NO

COMMENTS

None.

SIGNATURE

Clare Conley

2 of 8 cu

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-002265

07/05/2007

12:52:09

B4

Standard Lot#: 635301I

Last Changed By: T. GALLEGOS

#28847

QAS: C. CONLEY

#34517

TUCSON PD

Operator: C. CONLEY

#34517

TUCSON PD

Subject: TEST, SQAP

DOB: 11/11/1111

Sex: M

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	12:53:08
Diagnostic Test	Pass	12:53:36
Air Blank	0.000	12:54:03
0.100 Cal Check	0.100	12:54:22
Air Blank	0.000	12:54:51
Subject Test	0.000	12:55:12
Air Blank	0.000	12:55:39
Five Minute Wait	ABT*	
Air Blank	0.000	12:56:18

*Sequence Aborted

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-002265

07/05/2007

12:56:29

154

Standard Lot#: 635301I

Last Changed By: T. GALLEGOS

#28847

QAS: C. CONLEY

#34517

TUCSON PD

Operator: C. CONLEY

#34517

TUCSON PD

Subject: TEST, SQAP

DOB: 11/11/1111

Sex: M

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	12:57:28
Diagnostic Test	Pass	12:57:56
Air Blank	0.000	12:58:23
0.100 Cal Check	RFI*	12:58:34
Air Blank	0.000	12:59:04

*RFI Detect

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

4 of 8 Cu

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-002265

07/05/2007

12:59:28

TS4

Standard Lot#: 635301I

Last Changed By: T. GALLEGOS

#28847

QAS: C. CONLEY

#34517

TUCSON PD

Operator: C. CONLEY

#34517

TUCSON PD

Subject: TEST, SQAP

DOB: 11/11/1111

Sex: M

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	IPS*	13:00:14
Air Blank	0.000	13:00:42

*Improper Sample

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-002265

07/05/2007

13:01:12

Standard Lot#: 635301I

Last Changed By: T. GALLEGOS

#28847

QAS: C. CONLEY

#34517

TUCSON PD

Operator: C. CONLEY

#34517

TUCSON PD

Subject: TEST, SQAP

DOB: 11/11/1111

Sex: M

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	13:02:14
Diagnostic Test	Pass	13:02:42
Air Blank	0.000	13:03:09
0.100 Cal Check	0.100	13:03:28
Air Blank	0.000	13:03:57
Subject Test	0.082*	13:07:06
Air Blank	0.000	13:07:38
Air Blank	0.000	13:08:05
Subject Test	XXX**	13:08:24
Air Blank	0.000	13:08:58

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

b of 8 cu

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-002265

07/05/2007

13:09:50

Standard Lot#: 635301I


Last Changed By: T. GALLEGOS #28847

DH

QAS: C. CONLEY #34517

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	13:10:17
Cal Check	0.099	13:10:36
Air Blank	0.000	13:11:06

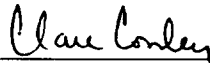

Operator's Signature

7/8/07

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

B4

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Clare Conley AGENCY Tucson Police Department
DATE 7-5-07 TIME 1250
Intoxilyzer Serial # 80-00 2265 Location 0DD

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 101 AC.
- () 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Clare Conley

DPS Form Exh G-2 (Rev. 05-01)

8 of 8 cu

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-002265

07/05/2007

12:49:40

64


Standard Lot#: 635301I

Last Changed By: T. GALLEGOS #28847

QAS: C. CONLEY #34517

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	12:50:07
Cal Check	0.101	12:50:26
Air Blank	0.000	12:50:56

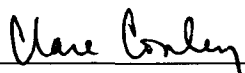

Operator's Signature