



CERTIFICATION by CUSTODIAN of RECORDS
Tucson Police Department Crime Laboratory

TS4
8-8-07

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002265 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 7-31-07 Cal 1 0.100 AC Cal 2 0.099 AC

Between the dates of 7-5-07 and 7-31-07 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.
The instrument was returned to service on N/A.


Clare Conley

Clare L. Conley

8-3-07

ADPS QAS Permit # 33314

SUBSCRIBED AND SWORN to before me this 3, day of August, 2007.

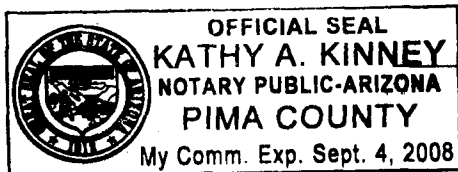

Kathy A Kinney
Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000B4
8-8-07

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Clare Conley AGENCY Tucson Police Department
DATE 7-31-07 TIME 13:50
INTOXILYZER SERIAL # 80-002265 LOCATION 0DD

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
(✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
(✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
(✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
(✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 099 AC

Instrument is operating properly and accurately. YES ☒ NO ☐

COMMENTS

None.

SIGNATURE

Clare Conley

2 of 8 Cu

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-002265

07/31/2007

13:46:14

Dh

Standard Lot#: 635301I

Last Changed By: T. GALLEGOS

#28847

QAS: C. CONLEY

#34517

TUCSON PD

Operator: C. CONLEY

#34517

TUCSON PD

Subject: TEST, SQAP

DOB: 11/11/1111

Sex: M

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	13:47:07
Diagnostic Test	Pass	13:47:35
Air Blank	0.000	13:48:02
0.100 Cal Check	0.100	13:48:21
Air Blank	0.000	13:48:50
Subject Test	0.000	13:49:09
Air Blank	ABT*	13:49:16
Air Blank	0.000	13:49:44

*Sequence Aborted

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

3 of 8 cu

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-002265

07/31/2007

13:52:19

B4

Standard Lot#: 635301I

Last Changed By: T. GALLEGOS

#28847

QAS: C. CONLEY

#34517

TUCSON PD

Operator: C. CONLEY

#34517

TUCSON PD

Subject: TEST, SQAP

DOB: 11/11/1111

Sex: M

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	13:53:19
Diagnostic Test	Pass	13:53:46
Air Blank	RFI*	13:53:55
Air Blank	0.000	13:54:24

*RFI Detect

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

4 of 8 cu

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-002265
07/31/2007

13:55:19

15h

Standard Lot#: 635301I
Last Changed By: T. GALLEGOS #28847

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP

DOB: 11/11/1111

Sex: M

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	IPS*	13:56:04
Air Blank	0.000	13:56:33

*Improper Sample

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

568 cu

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-002265

07/31/2007

13:56:59

ph

Standard Lot#: 635301I

Last Changed By: T. GALLEGOS

#28847

QAS: C. CONLEY

#34517

TUCSON PD

Operator: C. CONLEY

#34517

TUCSON PD

Subject: TEST, SQAP

DOB: 11/11/1111

Sex: M

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	13:57:54
Diagnostic Test	Pass	13:58:22
Air Blank	0.000	13:58:49
0.100 Cal Check	0.100	13:59:08
Air Blank	0.000	13:59:37
Subject Test	0.048*	14:02:46
Air Blank	0.000	14:03:18
Air Blank	0.000	14:03:46
Subject Test	XXX**	14:04:06
Air Blank	0.000	14:04:40

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-002265

07/31/2007

14:05:34

Standard Lot#: 635301I

Last Changed By: T. GALLEGOS #28847

QAS: C. CONLEY #34517

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	14:06:01
Cal Check	0.099	14:06:20
Air Blank	0.000	14:06:50

Clare Conley

Operator's Signature

7 of 8 m

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

156

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Clare Conley AGENCY Tucson Police Department
DATE 7-31-07 TIME 13:50
Intoxilyzer Serial # 80-002265 Location ODD

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0.100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0.100 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Clare Conley

DPS Form Exh G-2 (Rev. 05-01)

8 of 8 cu

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-002265

07/31/2007

13:50:07

ph

Standard Lot#: 635301I

Last Changed By: T. GALLEGOS #28847

QAS: C. CONLEY #34517

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	13:50:34
Cal Check	0.100	13:50:53
Air Blank	0.000	13:51:22

Clare Conley
Operator's Signature