



# CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
270 S. Stone Avenue Tucson, AZ 85701-1917

TS4  
3-11-2010

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-002265 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 03/11/10 Cal 1 0.100 AC Cal 2 0.099 AC

Between the dates of 02/11/10 and 03/11/10 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.  
The instrument was returned to service on N/A.

Cheryl K. Straub


03/11/10

ADPS QAS Permit # 33408

SUBSCRIBED AND SWORN to before me this 11<sup>th</sup>, day of March, 20 10.

Notary Public

\*Includes documents specified in A.A.C. R13-10-104.A.5

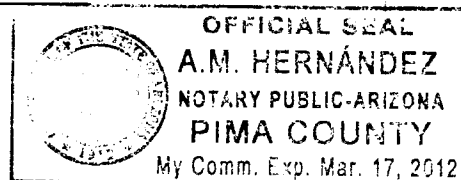

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

TBY  
3-11-2010QA SPECIALIST Cheryl STRAUB AGENCY Tucson Police DepartmentDATE 03/11/10 TIME 08:52INTOXILYZER SERIAL # 80-002265 LOCATION TEAM 5

(✓) 1. Display reads "PUSH BUTTON TO START"

## DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- 
- (✓) 2. Date check.

## OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0.
- 000
- AC.
- 
- (✓) 2. Error recognition logic system functioning.
- 
- Not a Successfully Completed Test Sequence printed
- 
- (✓) 3. Proper sample recognition system.
- 
- Not a Successfully Completed Test Sequence printed
- 
- Deficient sample printed.
- 
- (✓) 4. Standard Calibration Check standard 0.
- 100
- AC
- 
- Result: 0.
- 099
- AC

Instrument is operating properly and accurately. YES ✓ NO     

## COMMENTS

---

---

---

---

SIGNATURE

Cheryl Straub

297  
CS

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-002265

03/11/2010

08:53:52

BY

Standard Lot#: 642845

Last Changed By: C. STRAUB

#50772

QAS: C. STRAUB

#50772

TUCSON PD

Operator: C. STRAUB

#50772

TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	08:54:45
Diagnostic Test	Pass	08:55:13
Air Blank	0.000	08:55:40
0.100 Cal Check	0.100	08:55:59
Air Blank	0.000	08:56:29
Subject Test	0.000	08:56:45
Air Blank	ABT*	08:56:52
Air Blank	0.000	08:57:20

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

377  
CS

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-002265  
03/11/2010

08:57:24

754

---

Standard Lot#: 642845  
Last Changed By: C. STRAUB #50772

---

QAS: C. STRAUB #50772  
TUCSON PD

---

Operator: C. STRAUB #50772  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	RFI*	08:57:55
Air Blank	IPS**	08:58:03

\*RFI Detect  
\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

487  
CS

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-002265

03/11/2010

08:58:11

54

Standard Lot#: 642845

Last Changed By: C. STRAUB

#50772

QAS: C. STRAUB

#50772

TUCSON PD

Operator: C. STRAUB

#50772

TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	08:59:01
Diagnostic Test	Pass	08:59:28
Air Blank	0.000	08:59:55
0.100 Cal Check	0.100	09:00:14
Air Blank	0.000	09:00:43
Subject Test	0.102*	09:03:52
Air Blank	0.000	09:04:23
Air Blank	0.000	09:04:51
Subject Test	XXX**	09:05:15
Air Blank	0.000	09:05:50

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**

5 of 7  
CS

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-002265

03/11/2010

09:06:08

134

-----  
Standard Lot#: 642845

Last Changed By: C. STRAUB #50772  
-----

QAS: C. STRAUB #50772

TUCSON, PD  
-----

Test	g/210L	Time
Air Blank	0.000	09:06:35
Cal Check	0.099	09:06:54
Air Blank	0.000	09:07:23

  
-----

Operator's Signature

697  
CSTHIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000B4  
3-11-2010

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Cheryl STRAUB AGENCY Tucson Police Department  
DATE 03/11/10 TIME 08:52  
Intoxilyzer Serial # 80-002265 Location TEAM 5

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution N/A AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach  $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 100 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Cheryl Straub

DPS Form Exh G-2 (Rev. 05-01)

777  
CS

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-002265

03/11/2010

08:52:29

154


-----  
Standard Lot#: 642845

Last Changed By: C. STRAUB #50772  
-----

QAS: C. STRAUB #50772

TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	08:52:57
Cal Check	0.100	08:53:15
Air Blank	0.000	08:53:45

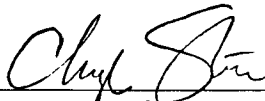

Operator's Signature