



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
270 S. Stone Avenue Tucson, AZ 85701-1917

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002265 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

TS4
5.6.2010

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 05/06/10 Cal 1 0.101 AC Cal 2 0.100 AC

Between the dates of 04/08/10 and 05/06/10 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.
The instrument was returned to service on N/A.

Cheryl K. Straub

Cheryl K. Straub

05/06/10


ADPS QAS Permit # 33408

SUBSCRIBED AND SWORN to before me this 6th, day of May, 20 10.

A.M. Hernández

Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

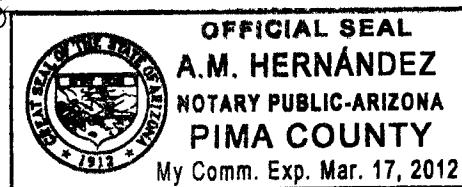

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

TS4
5.6.2010QA SPECIALIST Cheryl STRAUB AGENCY Tucson Police DepartmentDATE 05/06/10 TIME 09:38INTOXILYZER SERIAL # 80-002265 LOCATION TEAM 5

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

(✓) 1. Clock time check.

(✓) 2. Date check.

OPERATIONAL TESTS

(✓) 1. Alcohol-free subject test result 0. 000 AC.

(✓) 2. Error recognition logic system functioning.

Not a Successfully Completed Test Sequence printed

(✓) 3. Proper sample recognition system.

Not a Successfully Completed Test Sequence printed

Deficient sample printed.

(✓) 4. Standard Calibration Check standard 0. 100 ACResult: 0. 100 ACInstrument is operating properly and accurately. YES ✓ NO

COMMENTS

SIGNATURE

Cheryl Straub

277
CS

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-002265
05/06/2010

09:39:45

BH

Standard Lot#: 642845
Last Changed By: C. STRAUB #50772

QAS: C. STRAUB #50772
TUCSON PD

Operator: C. STRAUB #50772
TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:40:38
Diagnostic Test	Pass	09:41:06
Air Blank	0.000	09:41:33
0.100 Cal Check	0.101	09:41:52
Air Blank	0.000	09:42:21
Subject Test	0.000	09:42:47
Air Blank	ABT*	09:42:54
Air Blank	0.000	09:43:22

*Sequence Aborted

Not a Successfully Completed Test Sequence

377
CS

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-002265

05/06/2010

09:43:26

Standard Lot#: 642845

Last Changed By: C. STRAUB

#50772

QAS: C. STRAUB

#50772

TUCSON PD

Operator: C. STRAUB

#50772

TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	09:44:00
Air Blank	IPS**	09:44:09

*RFI Detect

**Improper Sample

Not a Successfully Completed Test Sequence

477
CS

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-002265

05/06/2010

09:44:16

134

Standard Lot#: 642845

Last Changed By: C. STRAUB

#50772

QAS: C. STRAUB

#50772

TUCSON PD

Operator: C. STRAUB

#50772

TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:45:08
Diagnostic Test	Pass	09:45:35
Air Blank	0.000	09:46:03
0.100 Cal Check	0.101	09:46:21
Air Blank	0.000	09:46:51
Subject Test	0.124*	09:49:59
Air Blank	0.000	09:50:31
Air Blank	0.000	09:50:59
Subject Test	XXX**	09:51:15
Air Blank	0.000	09:51:48

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

587
CS

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-002265

05/06/2010

09:52:06

TS4

Standard Lot#: 642845

Last Changed By: C. STRAUB #50772

QAS: C. STRAUB #50772

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	09:52:33
Cal Check	0.100	09:52:52
Air Blank	0.000	09:53:21



Operator's Signature

687
CSTHIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000B4
5.6.2010

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Cheryl STRAUB AGENCY Tucson Police Department
DATE 05/06/10 TIME 09:38
Intoxilyzer Serial # 80-002265 Location TEAM 5

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution N/A AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 101 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Cheryl Straub

DPS Form Exh G-2 (Rev. 05-01)

787
CS

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-002265

05/06/2010

09:38:24

64


Standard Lot#: 642845

Last Changed By: C. STRAUB #50772

QAS: C. STRAUB #50772

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	09:38:51
Cal Check	0.101	09:39:10
Air Blank	0.000	09:39:40

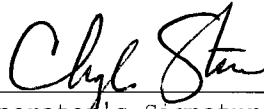

Operator's Signature