



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
270 S. Stone Avenue Tucson, AZ 85701-1917

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-002265 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

TSY  
7-6-2010

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 07/01/10 Cal 1 0.100 AC Cal 2 0.099 AC

Between the dates of 06/03/10 and 07/01/10 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.  
The instrument was returned to service on N/A.

Cheryl K. Straub

Cheryl K. Straub

07/01/10


ADPS QAS Permit # 33408

SUBSCRIBED AND SWORN to before me this 1, day of

July, 2010.

A.M. Hernandez

Notary Public

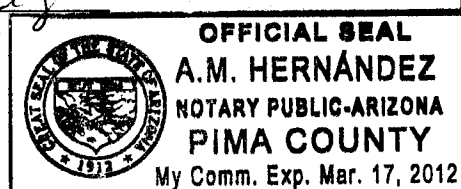

\*Includes documents specified in A.A.C. R13-10-104.A.5

187  
CS

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

By  
7-6-2010

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Cheryl STRAUB AGENCY Tucson Police Department

DATE 07/01/10 TIME 10:04

INTOXILYZER SERIAL # 80-002265 LOCATION TEAM 5

( ☒ ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- ( ☒ ) 1. Clock time check.  
( ☒ ) 2. Date check.

OPERATIONAL TESTS

- ( ☒ ) 1. Alcohol-free subject test result 0. 000 AC.  
( ☒ ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed  
( ☒ ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.  
( ☒ ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 099 AC

Instrument is operating properly and accurately. YES ☒ NO ☐

COMMENTS

\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE

Cheryl Straub

277  
CS

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-002265

07/01/2010

10:05:48

TS4

Standard Lot#: 642845

Last Changed By: C. STRAUB

#50772

QAS: C. STRAUB

#50772

TUCSON PD

Operator: C. STRAUB

#50772

TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

| Test            | g/210L | Time     |
|-----------------|--------|----------|
| Air Blank       | 0.000  | 10:06:39 |
| Diagnostic Test | Pass   | 10:07:07 |
| Air Blank       | 0.000  | 10:07:34 |
| 0.100 Cal Check | 0.100  | 10:07:53 |
| Air Blank       | 0.000  | 10:08:22 |
| Subject Test    | 0.000  | 10:08:39 |
| Air Blank       | ABT*   | 10:08:45 |
| Air Blank       | 0.000  | 10:09:14 |

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

317  
CS

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-002265  
07/01/2010

10:09:19

TS4

---

Standard Lot#: 642845  
Last Changed By: C. STRAUB #50772

---

QAS: C. STRAUB #50772  
TUCSON PD

---

Operator: C. STRAUB #50772  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

| Test      | g/210L | Time     |
|-----------|--------|----------|
| Air Blank | RFI*   | 10:09:56 |
| Air Blank | IPS**  | 10:10:05 |

\*RFI Detect

\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

477  
CS

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-002265  
07/01/2010

10:10:27

134

Standard Lot#: 642845  
Last Changed By: C. STRAUB #50772

QAS: C. STRAUB #50772  
TUCSON PD

Operator: C. STRAUB #50772  
TUCSON PD

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

| Test            | g/210L | Time     |
|-----------------|--------|----------|
| Air Blank       | 0.000  | 10:11:14 |
| Diagnostic Test | Pass   | 10:11:41 |
| Air Blank       | 0.000  | 10:12:08 |
| 0.100 Cal Check | 0.100  | 10:12:27 |
| Air Blank       | 0.000  | 10:12:57 |
| Subject Test    | 0.291* | 10:16:05 |
| Air Blank       | 0.000  | 10:16:39 |
| Air Blank       | 0.000  | 10:17:06 |
| Subject Test    | XXX**  | 10:17:37 |
| Air Blank       | 0.000  | 10:18:11 |

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**

577  
CS

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-002265

07/01/2010

10:18:31

184


-----  
Standard Lot#: 642845

Last Changed By: C. STRAUB #50772

-----  
QAS: C. STRAUB #50772

TUCSON PD  
-----

| Test      | g/210L | Time     |
|-----------|--------|----------|
| Air Blank | 0.000  | 10:18:59 |
| Cal Check | 0.099  | 10:19:17 |
| Air Blank | 0.000  | 10:19:47 |

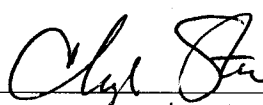  
-----

Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY  
STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Cheryl Straub AGENCY Tucson Police Department  
DATE 07/01/10 TIME 10:04  
Intoxilyzer Serial # 80-002265 Location TEAM 5

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution N/A AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach  $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 100 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Cheryl Straub

DPS Form Exh G-2 (Rev. 05-01)

777  
CS

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-002265

07/01/2010

10:04:24

154


Standard Lot#: 642845

Last Changed By: C. STRAUB #50772

QAS: C. STRAUB #50772

TUCSON PD

| Test      | g/210L | Time     |
|-----------|--------|----------|
| Air Blank | 0.000  | 10:04:51 |
| Cal Check | 0.100  | 10:05:10 |
| Air Blank | 0.000  | 10:05:39 |

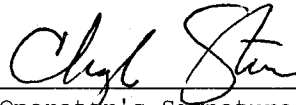

Operator's Signature