



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
270 S. Stone Avenue Tucson, AZ 85701-1917

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002245 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

by
8.4.2010

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 07/29/10 Cal 1 0.099 AC Cal 2 0.098 AC

Between the dates of 07/01/10 and 07/29/10 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.
The instrument was returned to service on N/A.

Cheryl K. Straub

Cheryl K. Straub

08/04/10

ADPS QAS Permit # 33408


SUBSCRIBED AND SWORN to before me this 4, day of

August, 2010.

A.M. Hernandez

Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

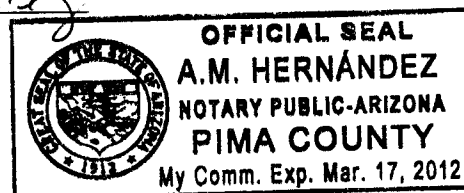

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

TS4
8.4.2010

QA SPECIALIST Cheryl STRAUB AGENCY Tucson Police Department
DATE 07/29/10 TIME 12:46
INTOXILYZER SERIAL # 80-002265 LOCATION TEAM 5

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

(✓) 1. Clock time check.
(✓) 2. Date check.

OPERATIONAL TESTS

(✓) 1. Alcohol-free subject test result 0. 000 AC.
(✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
(✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
(✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 099 AC

Instrument is operating properly and accurately. YES ✓ NO

COMMENTS

SIGNATURE

Cheryl Straub

DPS Form Exh G-4 (Rev. 05-01)

237
CS

Sh

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-002265
07/29/2010

12:48:07

Standard Lot#: 642845
Last Changed By: C. STRAUB #50772

QAS: C. STRAUB #50772
TUCSON PD

Operator: C. STRAUB #50772
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	12:48:58
Diagnostic Test	Pass	12:49:25
Air Blank	0.000	12:49:53
0.100 Cal Check	0.099	12:50:11
Air Blank	0.000	12:50:41
Subject Test	0.000	12:51:03
Air Blank	ABT*	12:51:09
Air Blank	0.000	12:51:38

*Sequence Aborted

Not a Successfully Completed Test Sequence

377
CS

B4

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-002265

07/29/2010

12:51:42

Standard Lot#: 642845

Last Changed By: C. STRAUB

#50772

QAS: C. STRAUB

#50772

TUCSON PD

Operator: C. STRAUB

#50772

TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	12:52:16
Air Blank	IPS**	12:52:27

*RFI Detect

**Improper Sample

Not a Successfully Completed Test Sequence

477
CS

154

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-002265

07/29/2010

12:52:35

Standard Lot#: 642845

Last Changed By: C. STRAUB

#50772

QAS: C. STRAUB

#50772

TUCSON PD

Operator: C. STRAUB

#50772

TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	12:53:27
Diagnostic Test	Pass	12:53:55
Air Blank	0.000	12:54:22
0.100 Cal Check	0.100	12:54:41
Air Blank	0.000	12:55:10
Subject Test	0.115*	12:58:19
Air Blank	0.000	12:58:51
Air Blank	0.000	12:59:19
Subject Test	XXX**	12:59:37
Air Blank	0.000	13:00:12

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

577
CS

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-002265

07/29/2010

13:00:30

TS4


Standard Lot#: 642845

Last Changed By: C. STRAUB #50772

QAS: C. STRAUB #50772

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	13:00:58
Cal Check	0.098	13:01:16
Air Blank	0.000	13:01:46

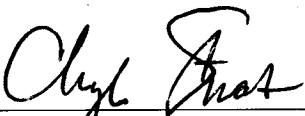

Operator's Signature

687
CS

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000


Dy
8.4.2010

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Cheryl STRAUB AGENCY Tucson Police Department
DATE 07/29/10 TIME 12:46
Intoxilyzer Serial # 80-002265 Location TEAM 5

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution N/A AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 099 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

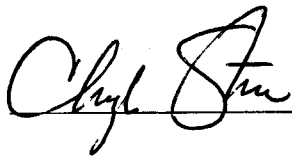

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-002265

07/29/2010

12:46:08


Standard Lot#: 642845

Last Changed By: C. STRAUB #50772

QAS: C. STRAUB #50772

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	12:46:36
Cal Check	0.099	12:46:54
Air Blank	0.000	12:47:24

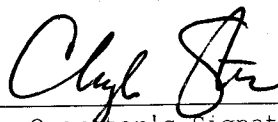

Operator's Signature