



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

TS4
9.27.2011

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002265 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 9/22/11 Cal 1 0.099 AC Cal 2 0.098 AC

Between the dates of 8/25/11 and 9/22/11 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.
The instrument was returned to service on -.

Sarah Cione

Sarah L. Cione


9/22/11

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 22, day of September, 20 11.

A.M. Hernandez

Notary Public

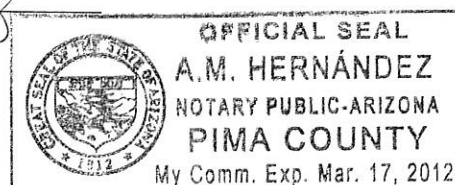

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

1 of 7sc

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

TJH
9.27.2011

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 9/22/11 TIME 1343

INTOXILYZER SERIAL # 80-002265 LOCATION Crime Lab

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
(✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
(✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
(✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
(✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 098 AC

Instrument is operating properly and accurately. YES ✓ NO

COMMENTS

SIGNATURE Sarah Cione

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-002265
09/22/2011

13:45:30

2 of 7sc

BH

Standard Lot#: 20910100A1
Last Changed By: C. STRAUB #50772

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	13:46:33
Diagnostic Test	Pass	13:47:01
Air Blank	0.000	13:47:28
0.100 Cal Check	0.099	13:47:47
Air Blank	0.000	13:48:16
Subject Test	0.000	13:48:34
Air Blank	0.000	13:49:01
Five Minute Wait	ABT*	
Air Blank	0.000	13:49:51

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-002265
09/22/2011

13:50:02

3 of 7sc
TS4

Standard Lot#: 20910100A1
Last Changed By: C. STRAUB #50772

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	13:50:51
Diagnostic Test	Pass	13:51:19
Air Blank	0.000	13:51:46
0.100 Cal Check	RFI*	13:52:01
Air Blank	IPS**	13:52:10

*RFI Detect

**Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-002265
09/22/2011

13:52:18

Standard Lot#: 20910100A1
Last Changed By: C. STRAUB #50772

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP,
DOB: 11/11/1111

Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	13:53:29
Diagnostic Test	Pass	13:53:56
Air Blank	0.000	13:54:23
0.100 Cal Check	0.099	13:54:42
Air Blank	0.000	13:55:12
Subject Test	0.018*	13:58:20
Air Blank	0.000	13:58:52
Air Blank	0.000	13:59:19
Subject Test	XXX**	13:59:45
Air Blank	0.000	14:00:18

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-002265

09/22/2011

14:00:47

Standard Lot#: 20910100A1

Last Changed By: C. STRAUB #50772

QAS: S. CIONE #50616

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	14:01:14
Cal Check	0.098	14:01:33
Air Blank	0.000	14:02:03

Sarah Cione

Operator's Signature

5 of 7 sc

BS

EXHIBIT G-2

6 of 7_{sc}

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

TSY
9-27-2011

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 9/22/11 TIME 1343
Intoxilyzer Serial # 80-002265 Location Crime Lab

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 099 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-002265

09/22/2011

13:43:59

Standard Lot#: 20910100A1

Last Changed By: C. STRAUB #50772

QAS: S. CIONE #50616

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	13:44:26
Cal Check	0.099	13:44:45
Air Blank	0.000	13:45:14

Sarah Cione

Operator's Signature