



# CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
1306 W. Miracle Mile, Tucson, AZ 85705

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-002265 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

TS4  
11-7-2011

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 10/19/11 Cal 1 0.100 AC Cal 2 0.099 AC

Between the dates of 9/22/11 and 10/19/11 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.  
The instrument was returned to service on -.

Sarah Cione

Sarah L. Cione


10/27/11

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 27, day of October, 2011.

A.M. Hernandez

Notary Public

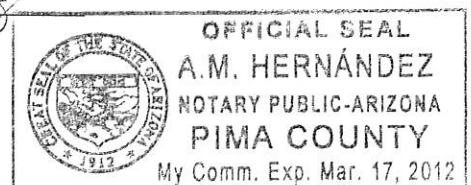

\*Includes documents specified in A.A.C. R13-10-104.A.5

1 of 7<sub>sc</sub>

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

TS4  
11-7-2011

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 10/19/11 TIME 1546

INTOXILYZER SERIAL # 80-002265 LOCATION Crime Lab

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.  
(✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.  
(✓) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed  
(✓) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.  
(✓) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 099 AC

Instrument is operating properly and accurately. YES ✓ NO     

COMMENTS

\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE Sarah Cione

2 of 7sc

134

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-002265

10/19/2011

15:48:44

Standard Lot#: 20910100A1

Last Changed By: C. STRAUB #50772

QAS: S. CIONE #50616  
TUCSON PDOperator: S. CIONE #50616  
TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	15:49:31
Diagnostic Test	Pass	15:49:58
Air Blank	0.000	15:50:26
0.100 Cal Check	0.100	15:50:44
Air Blank	0.000	15:51:14
Subject Test	0.000	15:51:31
Air Blank	0.000	15:51:58
Five Minute Wait	ABT*	
Air Blank	0.000	15:52:44

\*Sequence Aborted

### Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

3 of 7<sub>sc</sub>

TS4

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-002265

10/19/2011

15:52:53

---

Standard Lot#: 20910100A1

Last Changed By: C. STRAUB #50772

---

QAS: S. CIONE #50616

TUCSON PD

---

Operator: S. CIONE #50616

TUCSON PD

---

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	RFI*	15:53:29
Air Blank	IPS**	15:53:37

\*RFI Detect

\*\*Improper Sample

### Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

4 of 7<sub>sc</sub>  
T54

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-002265  
10/19/2011

15:54:45

---

Standard Lot#: 20910100A1  
Last Changed By: C. STRAUB #50772

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	15:56:04
Diagnostic Test	Pass	15:56:32
Air Blank	0.000	15:56:59
0.100 Cal Check	0.100	15:57:18
Air Blank	0.000	15:57:47
Subject Test	0.061*	16:00:56
Air Blank	0.000	16:01:28
Air Blank	0.000	16:01:55
Subject Test	XXX**	16:02:18
Air Blank	0.000	16:02:51

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

### Not a Successfully Completed Test Sequence

Pursuant to section 26-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

5 of 7<sub>sc</sub>

TS4

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-002265

10/19/2011

16:03:06

-----  
Standard Lot#: 20910100A1

Last Changed By: C. STRAUB #50772  
-----

QAS: S. CIONE #50616

TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	16:03:34
Cal Check	0.099	16:03:53
Air Blank	0.000	16:04:22

*Sarah Cione*

\_\_\_\_\_  
Operator's Signature

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

TS4  
11.7.2011

QA Specialist Sarah Cione AGENCY Tucson Police Department  
DATE 10/19/11 TIME 1546  
Intoxilyzer Serial # 80-002265 Location Crime Lab

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 100 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-002265

10/19/2011

15:46:58

7 of 7sc

Bl

-----  
Standard Lot#: 20910100A1

Last Changed By: C. STRAUB #50772  
-----

QAS: S. CIONE #50616

TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	15:47:25
Cal Check	0.100	15:47:44
Air Blank	0.000	15:48:14

Sarah Cione

Operator's Signature