



# CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
1306 W. Miracle Mile, Tucson, AZ 85705

TS4  
4.11.2012

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-002265 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 3/29/12 Cal 1 0.101 AC Cal 2 0.100 AC

Between the dates of 3/1/12 and 3/29/12 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.  
The instrument was returned to service on -.

Sarah Cione

Sarah L. Cione

4/9/12


ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 9th, day of April, 2012.

Elma Escarcega

Notary Public

\*Includes documents specified in A.A.C. R13-10-104.A.5

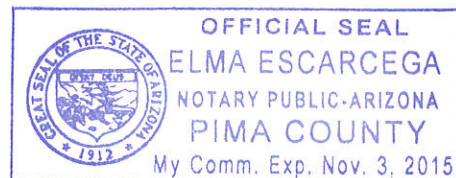

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

TS4  
4-11-2012

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 3/29/12 TIME 1343

INTOXILYZER SERIAL # 80-002265 LOCATION DUI HQ

( ☒ ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- ( ☒ ) 1. Clock time check.  
( ☒ ) 2. Date check.

OPERATIONAL TESTS

- ( ☒ ) 1. Alcohol-free subject test result 0. 000 AC.  
( ☒ ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed  
( ☒ ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.  
( ☒ ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 100 AC

Instrument is operating properly and accurately. YES ☒ NO ☐

COMMENTS

---

---

---

---

SIGNATURE Sarah Cione

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-002265

03/29/2012

13:45:50

Standard Lot#: 20910100A1

Last Changed By: C. STRAUB

#50772

QAS: S. CIONE

#50616

TUCSON PD

Operator: S. CIONE

#50616

TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	13:46:36
Diagnostic Test	Pass	13:47:04
Air Blank	0.000	13:47:31
0.100 Cal Check	0.101	13:47:50
Air Blank	0.000	13:48:20
Subject Test	0.000	13:49:02
Air Blank	0.000	13:49:29
Five Minute Wait	ABT*	
Air Blank	0.000	13:50:36

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

3 of 9 sc

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-002265  
03/29/2012

13:52:26

TS4

Standard Lot#: 20910100A1  
Last Changed By: C. STRAUB #50772

QAS: S. CIONE #50616  
TUCSON PD

Operator: S. CIONE #50616  
TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	13:52:55
Air Blank	IPS**	13:53:02

\*RFI Detect

\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

4 of 9<sub>sc</sub>

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-002265  
03/29/2012

14:00:09

B4

Standard Lot#: 20910100A1  
Last Changed By: C. STRAUB #50772

QAS: S. CIONE #50616  
TUCSON PD

Operator: S. CIONE #50616  
TUCSON PD

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	14:00:54
Diagnostic Test	Pass	14:01:21
Air Blank	0.000	14:01:48
0.100 Cal Check	0.101	14:02:07
Air Blank	0.000	14:02:37
Subject Test	0.136*	14:05:46
Air Blank	0.000	14:06:18
Air Blank	0.000	14:06:46
Subject Test	XXX**	14:07:15
Air Blank	0.000	14:07:49

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**



5 of 9<sub>sc</sub>

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 60-002265

03/29/2012

14:08:50

BS4


-----  
Standard Lot#: 20910100A1

Last Changed By: C. STRAUB #50772  
-----

QAS: S. CIONE #50616

TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	14:09:17
Cal Check	0.100	14:09:36
Air Blank	0.000	14:10:06

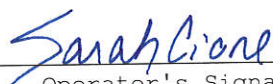

Operator's Signature

EXHIBIT G-2

6 of 9<sub>sc</sub>

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

TS4  
4.11.2012

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department  
DATE 3/29/12 TIME 1343  
Intoxilyzer Serial # 80-002265 Location DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 101 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-002265  
03/29/2012 13:43:27

7 of 9

TS4

-----  
Standard Lot#: 20910100A1  
Last Changed By: C. STRAUB #50772  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	13:43:55
Cal Check	0.101	13:44:14
Air Blank	0.000	13:44:44

  
\_\_\_\_\_  
Operator's Signature





TUCSON POLICE DEPARTMENT  
CRIME LABORATORY

8 of 95c

TS4

INTOXILYZER SERVICE RECORD  
FOR  
CALIBRATION GAS TANK CHANGE

4.11.2012

Date: 3/29/12  
QAS: Sarah Cione  
Intoxilyzer Serial Number: 80-002265  
New Lot Number: 04911100A1  
New Expiration Date: 3/1/2013

- (✓) 1. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 2. Remove the old gas tank and install the new one.
- (✓) 3. Hit the "ESC" key twice and swipe the QAS card.
- (✓) 4. Go to the "M-Maintenance" menu and select "S-Standard Change."  
Select "D" for dry gas and enter the following information:  
Enter Target Value: 0.100  
Enter Control Lot #: From the tank label  
Enter Exp Date: From the tank label (mm/dd/yyyy)  
Edit info Y/N? Y or N  
Update Std Info? Y ENTER  
"Please wait saving..."
- (✓) 5. Hit the "ESC" key twice.
- (✓) 6. Go to the "M-Maintenance" menu and select "T-Control Testing."  
Select "D" for dry control test. After selection is made press  
ENTER.
- (✓) 7. Air blank completed.
- (✓) 8. Calibration check completed. Test results 0. 100 AC.
- (✓) 9. Air blank completed.
- (✓) 10. Remove and sign printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-002265

03/29/2012

14:14:02


-----  
Standard Lot#: 04911100A1

Last Changed By: S. CIONE #50616  
-----

QAS: S. CIONE #50616

TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	14:14:30
Cal Check	0.100	14:14:48
Air Blank	0.000	14:15:19

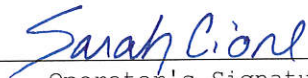

Operator's Signature