



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002265 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

TS4
6-27-2012

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 6/21/12 Cal 1 0.100 AC Cal 2 0.100 AC

Between the dates of 5/24/12 and 6/21/12 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.
The instrument was returned to service on -.

Sarah Cione

Sarah L. Cione

6/21/12


ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 21th, day of June, 20 12.

Edna Escarcega

Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

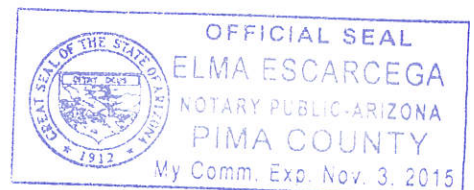

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

1 of 7_{sc}

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

TS4
627-2012

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department
DATE 6/21/12 TIME 1012
INTOXILYZER SERIAL # 80-002265 LOCATION DUT HQ

(☒) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (☒) 1. Clock time check.
(☒) 2. Date check.

OPERATIONAL TESTS

- (☒) 1. Alcohol-free subject test result 0. 000 AC.
(☒) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
(☒) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
(☒) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 100 AC

Instrument is operating properly and accurately. YES ☒ NO ☐

COMMENTS

SIGNATURE Sarah Cione

2 of 7_{sc}

TS4

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-002265

06/21/2012

10:14:22

Standard Lot#: 04911100A1

Last Changed By: S. CIONE

#50616

QAS: S. CIONE

#50616

TUCSON PD

Operator: S. CIONE

#50616

TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:15:03
Diagnostic Test	Pass	10:15:30
Air Blank	0.000	10:15:58
0.100 Cal Check	0.100	10:16:17
Air Blank	0.000	10:16:47
Subject Test	0.000	10:17:04
Air Blank	0.000	10:17:31
Five Minute Wait	ABT*	
Air Blank	0.000	10:18:14

*Sequence Aborted

Not a Successfully Completed Test Sequence

3 of 7sc

TS4

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-002265

06/21/2012

10:22:06

Standard Lot#: 04911100A1

Last Changed By: S. CIONE

#50616

QAS: S. CIONE

#50616

TUCSON PD

Operator: S. CIONE

#50616

TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	10:22:33
Air Blank	IPS**	10:22:42

*RFI Detect

**Improper Sample

Not a Successfully Completed Test Sequence

4 of 7_{sc}

Bh

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-002265

06/21/2012

10:25:42

Standard Lot#: 04911100A1

Last Changed By: S. CIONE

#50616

QAS: S. CIONE

#50616

TUCSON PD

Operator: S. CIONE

#50616

TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:26:31
Diagnostic Test	Pass	10:26:58
Air Blank	0.000	10:27:26
0.100 Cal Check	0.100	10:27:45
Air Blank	0.000	10:28:15
Subject Test	0.048*	10:31:23
Air Blank	0.000	10:31:55
Air Blank	0.000	10:32:22
Subject Test	XXX**	10:32:40
Air Blank	0.000	10:33:14

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-002265

06/21/2012

10:34:19

Standard Lot#: 04911100A1

Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	10:34:47
Cal Check	0.100	10:35:05
Air Blank	0.000	10:35:35

Sarah Cione

Operator's Signature

5 of 7 sc

Bh

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

6077_{sc}
TS4
6-27-2012

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 6/21/12 TIME 1012
Intoxilyzer Serial # 80-002265 Location DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 000 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 100 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

7 of 7_{sc}

BS4

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-002265

06/21/2012

10:12:17

Standard Lot#: 04911100A1

Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	10:12:44
Cal Check	0.100	10:13:03
Air Blank	0.000	10:13:32



Operator's Signature