



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

TS4
1-28-2013

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002265 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 1/4/13 Cal 1 0.100 AC Cal 2 0.100 AC

Between the dates of 12/7/12 and 1/4/13 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.
The instrument was returned to service on -.

Sarah Cione

Sarah L. Cione

1/7/13


ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 7th, day of January, 20 13.

Elma Escarcega

Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

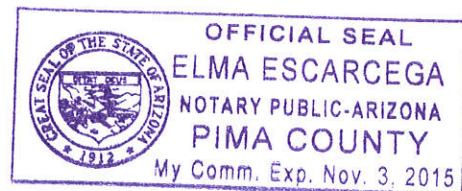

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1 of 7 sc
B4
1-28-2013

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 1/4/13 TIME 1453

INTOXILYZER SERIAL # 80-002265 LOCATION Crime Lab

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

(✓) 1. Clock time check.

(✓) 2. Date check.

OPERATIONAL TESTS

(✓) 1. Alcohol-free subject test result 0. 000 AC.

(✓) 2. Error recognition logic system functioning.

Not a Successfully Completed Test Sequence printed

(✓) 3. Proper sample recognition system.

Not a Successfully Completed Test Sequence printed

Deficient sample printed.

(✓) 4. Standard Calibration Check standard 0. 100 AC

Result: 0. 100 AC

Instrument is operating properly and accurately. YES ✓ NO

COMMENTS

SIGNATURE Sarah Cione

2 of 7sc
Dy

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-002265
01/04/2013 15:04:36

Standard Lot#: 09312100A1
Last Changed By: S. CIONE #50616

OAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SOAP, SOAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	15:05:19
Diagnostic Test	Pass	15:05:48
Air Blank	0.000	15:06:15
0.100 Cal Check	0.101	15:06:34
Air Blank	0.000	15:07:04
Subject Test	0.000	15:07:30
Air Blank	0.000	15:07:57
Five Minute Wait	ABT*	
Air Blank	0.000	15:08:42

*Sequence Aborted

Not a Successfully
Completed
Test Sequence

3 of 7_{sc}
TX

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-002265
01/04/2013 15:09:01

Standard Lot#: 09312100A1
Last Changed By: S. CIONE #50616

OAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SOAP, SOAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	15:09:44
Diagnostic Test	Pass	15:10:13
Air Blank	0.000	15:10:40
0.100 Cal Check	RFI*	15:10:55
Air Blank	IPS**	15:11:03

*RFI Detect
**Improper Sample

Not a Successfully
Completed
Test Sequence

4 of 7_{sc}
TS4

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-002265
01/04/2013 15:11:11

Standard Lot#: 09312100A1
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	15:11:55
Diagnostic Test	Pass	15:12:24
Air Blank	0.000	15:12:52
0.100 Cal Check	0.099	15:13:10
Air Blank	0.000	15:13:40
Subject Test	0.122*	15:16:49
Air Blank	0.000	15:17:21
Air Blank	0.000	15:17:48
Subject Test	XXX**	15:18:35
Air Blank	0.000	15:19:09

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully
Completed
Test Sequence

5 of 7_{sc}
TS4

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-002265
01/04/2013 15:21:14

Standard Lot#: 09312100A1
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	15:21:41
Cal Check	0.100	15:22:00
Air Blank	0.000	15:22:30

Sarah Cione
Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 1/4/13 TIME 1453
Intoxilyzer Serial # 80-002265 Location Crime Lab

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 100 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

6 of 7sc

BH
1.28.2013

7 of 7_{sc}
B4

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-002265
01/04/2013 14:53:22

Standard Lot#: 09312100A1
Last Changed By: S. CIONE #50616

OAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	14:53:50
Cal Check	0.100	14:54:08
Air Blank	0.000	14:54:38

Sarah Cione

Operator's Signature